

Protocol for Endoscopic Procedures – Dr. Shakir


For all procedures:

- I prefer the Ambu flexible scope
- **No urinalysis needed** unless the patient has symptoms of a UTI
- Antibiotics: Cipro x 1 tablet unless the patient has an allergy or is already on antibiotics

Procedure	Indication	Special equipment	Prep/Drape
Standard cystoscopy	• Hematuria	• 60cc syringe cytology	Genitalia
	• Evaluate prostate	• None	
	• URS stent pull	• Stent grasper • Tuohy-Borst/UroLock/SureSeal	
Cysto+urethral stricture dilation	Known or suspected urethral stricture	• Boston Scientific UroMax 24 French x 4cm balloon dilator(Boston Scientific Encore inflator* • Saline or water • Sensor 0.035 straight wire • Tuohy-Borst/UroLock/SureSeal	Genitalia
Neobladder Endoscopy	Stent pull after neobladder (recent surgery)	• Stent grasper • Tuohy-Borst/UroLock/SureSeal	Genitalia
	Evaluate neobladder (prior surgery)	None	
Ileal Conduit Endoscopy	Stent pull after ileal conduit (recent surgery)	• Stent grasper • Tuohy-Borst/UroLock/SureSeal	Stoma
	Evaluate ileal conduit (prior surgery)	None	

*Product number M0062251040. We should have at least 3 full kits (dilator + inflator) available in clinic at all times. Sometimes they may get used on call.

Setting Up UroMax Dilator

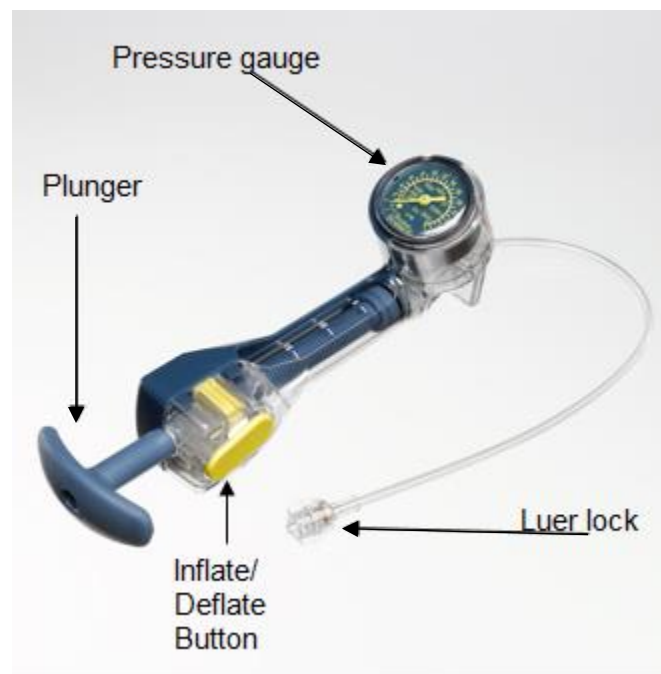

Setting Up UroMax Dilator (continued)

1. Open the balloon and inflator after confirming that we will be doing a urethral stricture dilation.
2. Use the box the inflator came in to pour some saline or water.
3. *Sterile*: Put the Luer lock end of the inflator into the saline. Push the yellow button, and **while keeping the button pushed**, pull back on the plunger. The inflator will fill with fluid.
4. *Sterile*: Connect the Luer lock end of the inflator to the “**Balloon**” end of the balloon dilator.
5. Cystoscopy is done with the Tuohy-Borst on the cystoscope. Once the stricture is identified, the 0.035 straight Sensor wire is placed through the scope and through the stricture into the bladder. The cystoscope is removed.
6. *Sterile*: The UroMax balloon is placed over the wire. The cystoscope is passed alongside it.
7. The balloon is inflated under cystoscopic vision:
 - a. **Setting the pressure**: Push the yellow button, and **while keeping the button pushed**, push the plunger in slightly. Then, **release the yellow button**.
 - b. Start turning the plunger clockwise. The needle on the pressure gauge will rise and the balloon should start inflating.
 - c. Inflate to 20 atm and hold for about 30 seconds.
 - d. To deflate, push the yellow button, **and while keeping the button pushed**, pull back all the way on the plunger until the fluid returns in the inflator, and the balloon is deflated.
 - e. Repeat treatments as necessary.

Video of the above:

<https://twitter.com/i/status/1330330380575105028>

